

PAB[®] Container

Helpful Hints

Refer to product labeling for full prescribing information.

Available in 25 mL fill in 100 mL (116 per case); 50 mL fill in 100 mL (84 per case); and 100 mL fill in 150 mL (64 per case). The containers are packaged on their side.



ADDING MEDICATION

Step 1: Remove the blue protective cap. *If the blue protective cap is missing, set aside and do not use.*


Step 2: Prior to adding medications, swab injection port. Allow to air dry.

Step 3: Support the base of the medication port, fully insert needle and inject medication through the center of injection site keeping fingers behind the flange. Remove needle. Mix thoroughly by inverting the container several times.

SPIKING

Step 4: Remove the foil seal on the set port in two steps: lift the tab from side of port; peel off towards additive port.

Step 5: Hold the flange of the set port between your thumb and forefinger. Push spike straight into port while twisting.

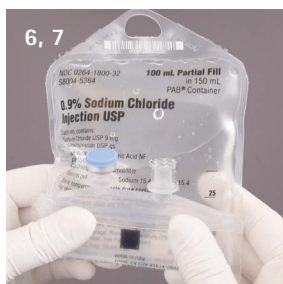

PNEUMATIC TUBE TRANSPORT

Step 6: Fold PAB IV container in half so that the hanger tab is positioned on top and extends past the ports to cushion the ports.

Step 7: Load PAB IV container into the pneumatic tube carrier in the folded half position with the hanger tab positioned as shown.

Step 8: Ensure the carrier is securely latched; watch for pinch points.

PNEUMATIC TUBE TRANSPORT WITH ADDEASE[®]

Step 9: Gently fold PAB diagonally; use caution not to place extra pressure on the container to avoid premature activation.

Step 10: Keeping the diagonal fold on the PAB, position the connection in the pneumatic tube carrier with the vial facing down. Close carrier carefully to ensure there are no pinched sides. Close latches.



SCANNING

Hold the center of the container parallel to your body and scan at a 45°-90° degree angle, 2"-6" from the container.

LOT NUMBER AND EXPIRATION DATE

The expiration date and lot numbers are embossed on the hanger tab at the top of the container.



B. Braun Medical Inc. | Bethlehem, PA
800-227-2862 | www.BBraunUSA.com